

Preservation above-zero

Blood Bank Refrigerators

Specially designed to provide precise and reliable temperature control that meets various criteria for whole blood and blood derivatives.



MBR-506D(H)

The Panasonic MBR Series provides two styles of storage systems: G-type (coated hard steel wire shelf) and GR-type (stainless steel roll-out drawers with card holder).

Designed to conform to AABB criteria, the series ensures stable and reliable temperature control utilizing Panasonic original technology. A special highly efficient compressor designed and developed by Panasonic provides rapid cooling and quiet performance for each model.

Temperature Variations Prevented

The Panasonic MBR Series is designed to minimize cold air loss even with frequent door openings.

- **Separated transparent inner doors** minimize the chamber air leakage during door openings.
- **Foamed-in-place insulation** in the walls and **magnetic sealed outer doors with double-pane glass window** prevent chamber air leakage and promote complete door closings.
- **Large air circulation fan** enables rapid temperature recovery after door openings.

User-Friendly Design

- **Selectable storage system** (G-type and GR-type) to suit user's needs.
- **Fluorescent interior lamp** with ON/OFF switch and a **large view window** in the outer door provide a clear view of stored items.
- **Digital display** is easy to see, and is calibratable through the control panel.



Easy-access front control panel (MBR-305GR)

- Filterless construction eliminates bothersome filter cleaning.

Alarm and Safety Functions


To ensure the safety of precious blood supply, the Panasonic MBR Series provides the following safety functions:

- **Audible and flashing LED visual alarms** with remote alarm contacts, in case of power failure, high or low temperature condition, or any thermal sensor abnormality.
- Door alarm and key lock are standard features.
- Re-activating buzzer, lamp and remote alarm contact. (30 min. after buzzer stops)
- Built-in temperature recorder.

Environmentally Friendly

Panasonic refrigerators feature commercially-available HCFC-free, HFC refrigerants and HCFC-free insulation.

Dimensions Unit : mm (inch)

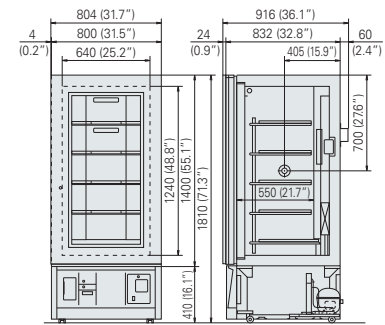

MBR-506D(H)

SPESIFICATION

Model No.	MBR-506D(H)
Exterior dimensions (W x D x H)	800 x 832 x 1810 (mm) / 31.5 x 32.8 x 71.3 (inch)
Interior dimensions (W x D x H)	640 x 550 x 1240 (mm) / 25.2 x 21.7 x 48.8 (inch)
Net weight	185kg
Effective capacity	425 liters [15.0 cu.ft.]
Storage capacity	120 bags (450ml), with MBR-56B
Exterior cabinet	Galvanised steel with baked-on finish
Interior cabinet	Stainless steel
Cabinet insulation	Foamed-in-place rigid polyurethane
Outer door	Triple layer glass windows
Door lock	1
Inner door (acrylic)	5
Shelves	5 coated hard steel wire shelves
Casters	4
Compressor	160W hermetically sealed
Refrigerant	HFC
Temperature	4°C ±1°C (Ambient temperature 35°C)
Air circulation	Forced air circulation
Defrost	Fully automatic
Temperature control	Microprocessor
Temperature alarm	High (6°C), Low (2°C), Audible and visual alarm
Power failure alarm	Audible and visual alarm (9 hours), Automatic rechargeable battery (Ni-Cd)
Door alarm	Visual alarm
Remote alarm contact	DC24V, 1A
	Temperature alarm or power failure alarm
Lighting	15W fluorescent lamp
Access port	1
Accessories	30-day recorder, Door lock key 1 set
Options	Recording paper (RP-06-PW)
	Basket (MBR-55B-PW) Max. 20pcs/unit 220ml bags x 10pcs/basket
	Basket (MBR-56B-PW) Max. 20pcs/unit 450ml bags x 6pcs/basket

Caution: Panasonic guarantees the product under certain warranty conditions. Panasonic in no way shall be responsible for any loss of content or damage to content.
• Appearance and specifications are subject to change without notice.



Panasonic Healthcare Co., Ltd., Gunma Factory is certified for:
Quality management system: ISO9001
Medical devices quality management system: ISO13485



Panasonic Healthcare Co., Ltd., Gunma Factory is certified for:
Environmental management system: ISO14001

DISTRIBUTED BY:

PT. SIGMA BIMED
Indonesia Sole Distributor

Jl. Agung Utara Raya Blok A 36D/42 Sunter Agung Jakarta Utara 14350 Indonesia

Telp. (62-21) 64716566, 70240088, 70240099, 87960123, 87960456

Fax. (62-21) 64716566, 87961123, 87961456

Email : sigmabimed@indo.net.id, sales@sigmabimed.com

Website : www.sigmabimed.com
kemenkes RI AKL

Panasonic

Panasonic Healthcare., Ltd.

1-1-1 Sakata, Oizumi-Machi, Ora-Gun, Gunma 370-0596, Japan
©Panasonic Healthcare Co., Ltd. 2012